

STATCLEAN™

Ultrasonic Cleaner



Code: **SC-P4862**

TREAT YOURSELF TO A QUIET ULTRASONIC

StatClean's 1.6 gallon tank and convenient self docking lid, holds a variety of instruments in the unique hands free drying basket, making it the perfect ultrasonic to complete your reprocessing needs.



Quiet

StatClean's 4.5 mm thick wall housing, made from engineering grade plastic, is quieter than most steel housing designs. This makes StatClean one of the quietest ultrasonic units on the market.



Easy to Use

StatClean's 3 color LED display accurately displays the 1 – 30 minute full range timer and Degas* option. The StatClean is equipped with advanced power-off memory, allowing it to save the last timer setting.

The easy to access drainage valve eliminates the need to manually empty a tank of solution, eliminating the possibility of spills.



Power

Two industrial grade ultrasonic transducers with independent control circuits provide uniform distribution of ultrasonic waves, increasing cleaning efficiency by up to 40% over other units on the market.



Water Resistant


The StatClean has been designed to be resistant to water penetration within the control panel to ensure continued optimal performance.

*Please consult your instrument cleaning solution instructions for use to determine if degassing is required. OPTIM 28 does not require degassing.

StatClean Specifications

Unit Size	20.3" (51.5 cm) x 12.0 (30.5 cm) x 9.4" (24.0 cm)
Power Supply	260W (AC 100-120V, 50/60 HZ)
Tank Capacity	1.6 gallon (6.0L)
Maximum Capacity	1.3 gallon (5.0L)
Minimum Capacity	0.9 gallon (3.35L)
Tank Size	12.7" x 8.1" x 3.9" (32.3 x 20.5 x 9.9 cm)
Longest item fits inside tank	13.4" (34.0 cm)


Accessories

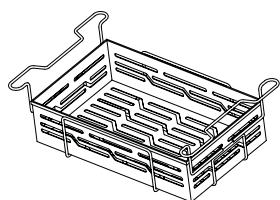

Metal Basket
SC-SB4862A



Power Cord
SC-PC4862A

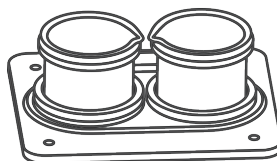

Double Beaker Holder
SC-BH02A



Hose
SC-PH01A



OPTIM[®] 28[†]

Ultrasonic Cleaning Solution

OPTIM[®] 28 is the cleaner of choice for use in the StatClean[™] and provides quick, superior cleaning of instruments. OPTIM[®] 28 contains Stabilized Hydrogen Peroxide, a highly effective means of cleaning organic debris. It is non-corrosive when used at the correct dilution and is effective against protein, carbohydrates, and lipids in only 5 minutes.

†Manufactured by:
Virox Technologies Inc.
2770 Coventry Rd, Oakville, ON
Canada, L6H 6R1
Tel: (905) 813-0110

Manufactured for + Distributed by:
SciCan Ltd.
1440 Don Mills Rd., Toronto, ON,
Canada, M3B 3P9
Tel: 416-445-1600
Toll free: 1-800-667-7733
customerservice@scican.com

Distributed by:
SciCan Inc.
701 Technology Dr., Canonsburg,
PA 15317, USA
Tel: 724-820-1600
Toll Free: 1-800-572-1211
uscustomerservice@scican.com